
AFM-2821


Aprilla

EN Instruction Manual
FR Mode D'emploi
DE Bedienungsanleitung
TR Kullanım Kılavuzu

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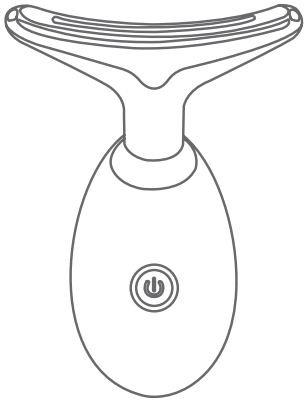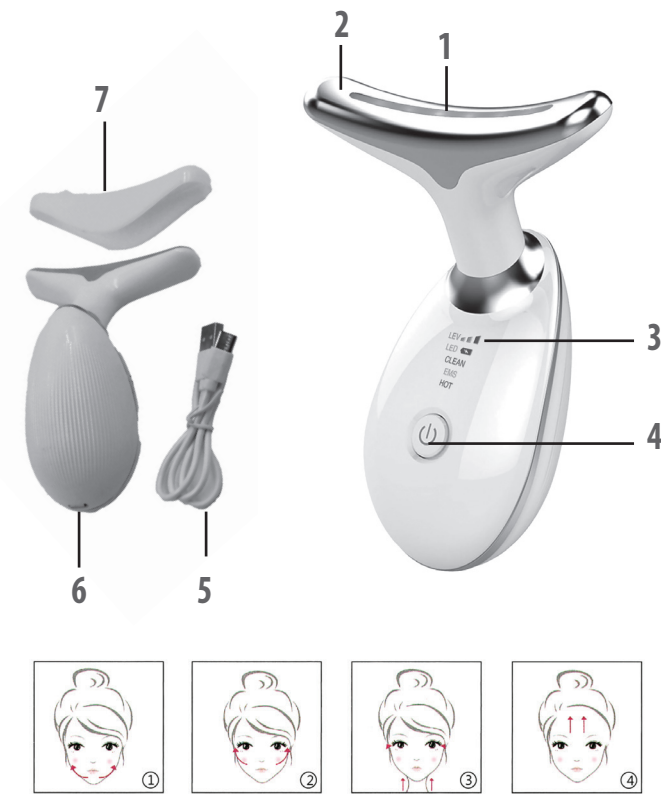EN Instruction Manual

SAFETY
PLEASE READ THE IMPORTANT INSTRUCTIONS FOR USER SAFETY CAREFULLY
AND KEEP FOR FUTURE REFERENCE!

- The manufacturer cannot be held liable for damage if the safety instructions are not followed.
- If the power cord is damaged, it must be replaced by the manufacturer, its service agent or a similarly qualified person in order to avoid a hazard.
- Never move the appliance by pulling on the cord and make sure that the cord cannot get tangled.
- The appliance must be placed on a stable, flat surface.
- The user must not leave the appliance unattended when connected to the power supply.
- This appliance must only be used in private households and only for the purpose for which it was designed.
- To protect yourself from electric shock, do not immerse the cord, plug or appliance in water or any other liquid.
- The appliance must only be connected to an earthed socket. Before connecting the plug to the socket, make sure that the voltage marked on the appliance corresponds to the local mains voltage.

CHILD SAFETY: This appliance must not be used by children under 8 years of age. This device may be used by children over 8 years of age and by people with reduced physical, sensory or mental capabilities or lack of experience and knowledge only if they are supervised or understand the hazards and instructions for safe use of the device. Children should not play with the device. Keep the device and its cord out of reach of children under 8 years of age. Cleaning and maintenance should not be performed by children unless they are over 8 years of age and under supervision.

GENERAL INFORMATION: EMS is the abbreviation of electromagnetic sound transduction and electrical muscle stimulation. It is a magnetic sound burst that strengthens the muscles in your face and increases the structure of the muscle. Radiofrequency is a heat method. It provides firmness to the skin by stimulating the collagen and fibers that allow the most common layer of the skin to stretch. With the distribution of heat, the collagen fibers begin to contract and renew, increasing the tightening and stretching of the skin.

WARNING: EMS devices should be used no more than 3 times a week and at 48-hour intervals.

DEVICE SAFETY:

- Do not immerse this electronic device, charging adapter or USB cable in water, do not wash in running water.
- Do not use or charge the product near running water or places containing water. Do not store the product in humid places such as bathrooms.
- If you drop this electronic device in water while charging, never touch it. Immediately unplug the charging adapter or USB cable from its socket. Make sure that all components are dry before operating it again.
- Do not touch the USB port, adapter or socket with wet hands while charging this electronic device.
- This product has LED (light emitting diode) feature. Never look directly at the light source with the naked eye or hold the device in your eyes.
- Do not use, charge or store this electronic device near a kitchen stove or in hot and dusty places such as a car. This may cause the product to malfunction or shorten its life.

Do not use this device in the following cases.

- If you are undergoing any skin treatment or have a skin condition.
- If you have extremely tanned skin.
- If you have bleeding under the skin that causes your skin to appear red or purple.
- If you have severely inflamed, open sores or itchy, red or rashy acne.
- If you have a light allergy or are taking a medication that can cause a light allergy (such as tetracycline).
- If you have a metal allergy.
- If you have a skin wound.
- If you have itchy, dry, rashy and irritated skin.
- If you are under 18 years of age.
- If you are receiving radiotherapy or chemotherapy.
- If you have a serious heart condition or use a pacemaker.
- If you have a thyroid disease or epilepsy.
- If you have a blood circulation disorder or diabetes, tuberculosis, lupus, a tumor in your body, a neurological disease that causes severe pain in your face (facial neuralgia) or a genetic disorder that causes hypersensitivity to light (porphyria).
- If you have a bone fracture or osteoporosis in your face.
- If you are taking a hormone medication.
- If you are pregnant, breastfeeding or menstruating.
- If you have previously had a skin treatment such as Botox or fillers applied to your skin, if you have had plastic surgery or if you are unsure about anything, consult your doctor before using this product.
- Do not use this product on your eyes, ears, lips, inside your mouth or private areas.

WARNING: The following situations may cause possible injuries and serious material damage, so you need to be very careful.

- Do not use the product around or near your eyes.
- Do not share the same product with anyone else.

- Do not use this product if you cannot perceive heat correctly. Diseases such as diabetes can disrupt the body's heat perception system. Disabled or sick people may not be able to express when they feel too hot. For this reason, do not use the product on these people.
- The device has a temperature protection system, but still, do not hold your device on the same part of your face for a long time (more than 3 minutes) while massaging, as this may cause burns.
- The recommended frequency of use may vary depending on the person's skin type or condition. Adjust your frequency of use according to your skin condition.
- If you feel any discomfort while using the product or after use, stop using it immediately and consult a doctor.

PARTS DESCRIPTION

- | | |
|---------------------------|-------------------------------|
| 1. Blue, red, green light | 2. Press sensitive head |
| 3. Speed indicator | 4. ON/off, mode select button |
| 5. Charge cable | 6. Charge input |
| 7. Protection cap | |

BEFORE FIRST USE

- Remove the device and its accessories from the box. Remove any labels, protective foil or plastic from the device.
- Before using your device for the first time, wipe all removable parts with a damp cloth. Never use abrasive products.
- Before first use, charge the device for at least 2 hours.
- To charge the device, connect the USB cable in the box to an adapter (or external battery) with a voltage higher than DC 5V or a current output of 1.5A. (Using an unsuitable adapter may cause the device to malfunction.)
- When you connect it to the USB cable, you can see the charging status from the charging indicator light on the device.

USING

1.Using Methods of Cleaning and exporting mode

- * This mode can effectively clean the dirt in the pores and make the skin become translucent; please pay attention to not using too thick cotton pads.
- * Combined with areen liiht. it can improve cells/olternative function of oxygen utilization, facilitate skin microcirculation, accelerate lymph circulation and reduce edema, and calm the skin to a stable condition.
- (1) Hold down button for about 1 second to enterthe "CLEAN" mode (accompanied with green light)
- (2) Apply makeup remover or other cleansing lotions onto the skin equably (instrument body is not waterproof, please don't directly spray water on it)
- (3) Hold the metal induction bars on the two sides (forming a loop with human body to realize the realion import and export)
- (4) Lift up the instrument head on the face or neck from bottom lo top

2. Using Methods of EMS nursing mode

- * This mode can effectively stimulate and relax the skin and help skin care products moisten the deep skin; cotton pad fixing ring can be removed. (There is a slight sense of current)
- *Combined with red light, it can effectively penetrate into deep skin, revitalize collagen, remove darkness and reduce fine lines and wrinkles.
- (1) Hold down function button for about 1 second to start and press function button again to enter the "EMS" mode (accompanied with red light)
- (2) Apply toner or lotion for daily use onto the skin equably
- (3) Hold the metalinduction bars on the two sides (forming a loop with human body to realize the real ion import and export)
- (4) Lift up the instrument head on the face or neck from bottom to top

3. Using Methods of Warm import mode

- * This mode can gently wake and relax skin and help skin recover lustre and creativity.
- * Combined with bule light, it can facilitate the collagen synthesis, rejuvenates skin, tighten sagging skin, refine skin pores and obvious improvements can be seen on the oily skin, sensitive skin and acne skin.
- (1) Hold down function button for about 1 second to start and press twice button to enter the "HOT" mode laccompanied with blue liiht)
- (2) Apply mask on the face or neck
- (3) Hold the metal induction bars on the two sides iformina a loop with human body to realize the rea! ion import and export)
- (4) Lift up the instrument head on the face or neck from bottom to top

Recommended usage (Side figures)

1. Lift up from chin to cheek
2. Lift up from cheek to outer canthus
3. Lift up from neck to chin
4. Lift up from chin to maxilla

CLEANING AND MAINTENANCE

- WARNING: Before performing any maintenance on the unit, make sure the device is turned off and unplugged.
- Clean the device with a damp cloth. Never use harsh and abrasive cleaners, scouring pads or steel wool, which may damage the device.
 - Never immerse the electric device in water or any other liquid.

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